


Daylight ID-Printer System

(part# ID-B2)

The Colenta ID Printer is for both AP and PA applications and incorporates the latest technology to record all the important patient data clearly and precisely onto each x-ray through the window area of the cassette in daylight operating conditions.




ID Card Dimensions

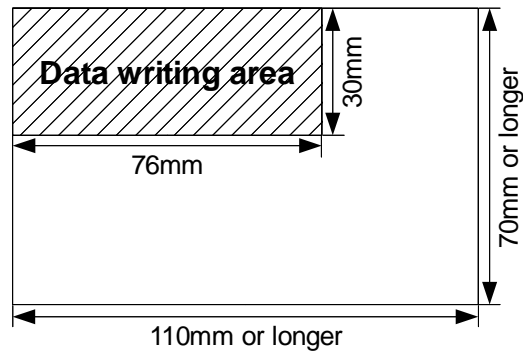

* Please use ID card with dimensions of 70 mm x 110 mm or larger to suit ID card slot

Type.....AP/PA Daylight ID printing system used with window cassette
 Functions.....prints patient's data written or typed on the white card which user provides.
 Month, date and year with time are automatically printed.
 Printing Method.....by tungsten bulb flashing
 Required Cassette.....window type cassette
 Density Adjusting.....density is adjustable in 4 steps, indicating by LED
 Safeguard.....The printer works when ID card is properly inserted into the card holder.
 4 types error indication
 Clock.....built-in digital clock is battery backup. It prints month, date, year and time.
 Power Source.....single phase, AC 100-240V, 50/60Hz
 Measurements and Weight...295 mm (L) x 420 mm (W) x 310 mm (H), about 10 kg

Printed sample on film

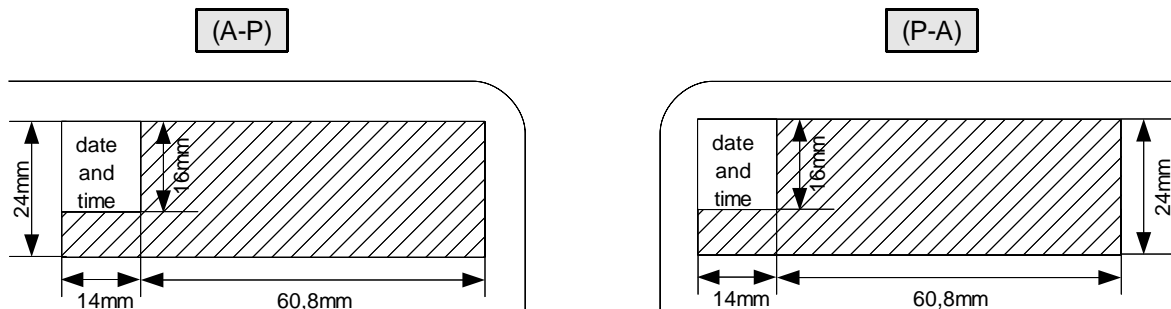

Printing dimensions are 24 mm x 60,8 mm - plus date and time.