

Colenta[®] MEDIPHOT 902/903

COLENTA MEDIPHOT X-Ray Film Processor: TYPE: MP 902 and MP 903

For processing all regular used type of Cut sheets of medical x-ray films.


	DEV-Time	Cycle Time	Intake Speed
MP 902 (max 43cm)	16sec to 90 sec	80 sec 440 sec	150 cm/min to (min.) 30 cm/min (max.)
MP 903 (max 43cm)	19sec to 105 sec	81 sec 450 sec	180 cm/min to (min.) 33 cm/min (max.)

All MEDIPHOT 902/903-models have 9 programmable processing channels, an intermediate rinse bath rack system, feed in table, hose/tube connection (100 mm dia.) facility for ventilation of the processors wet-section to an external exhaust fan and supplied with 2 x 30 lt replenishment storage tanks





Standard Features and accessories:

- Low Level monitoring of all Tank solutions.
- Automatic Solution Tank „Top-Up“ control on „Start-Up“
- „Work in Progress“ monitoring
- Rep storage tanks
- prepared to connect the processor to an external exhaust air system

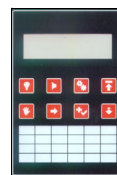

Microprocessor controlled Display

TRANSPORT SYSTEM:

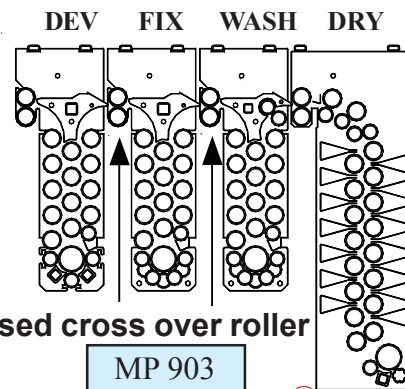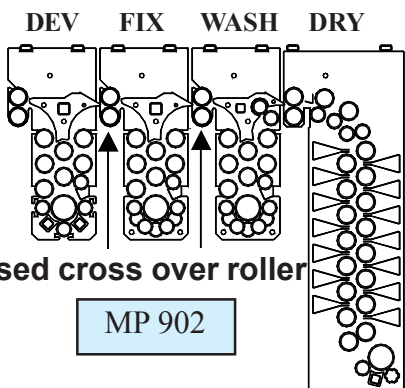

Pull-out tray for easy handling the replenishment containers




Optional accessories:

Cover Lid for Film Feed Tray



Rack Carrier Tray/Drip Tray for easy handling the transport racks

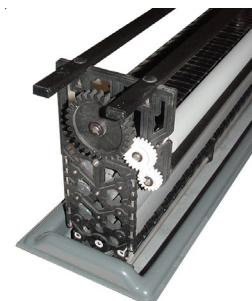

Colenta

Technical specification subject to change without notice.

see additional Datasheet X1-14

COLENTA LABORTECHNIK
Ges.m.b.H. & Co. KG
Neunkirchner Strasse 117
A-2700 Wiener Neustadt
Tel.: +43-2622-28311-0
Fax: +43-2622-28311-7
E-mail: office@colenta.at
<http://www.colenta.at>

ACI-COLENTA
INTERNATIONAL Div.
Forge Close, Little End Rd,
Eaton Socon, Cambs. UK
Tel. +44-1480-211118
Fax: +44-1480-215236
Fax: +44/870/132 1352
E-mail: aci@patflaherty.com