

MEDIPHOT 937/943 X-RAY



COLENTA Medical X-Ray Film Processors: MP937 and MP943

For all standard X-Ray Films - Certified for Mammography - Film processing.

	DEV-Time	Cycle Time	Intake Speed
MP937/943:	18sec to	1'24"min to	106cm/min to (min.)
max. 37 (43) cm	99 sec	7'43"min	20cm/min (max.)

All MEDIPHOT-model processors have 9 programmable processing channels, an intermediate rinse bath rack system, film feed tray (optional with Cover lid) and supplied with 2 x 30 lt replenishment storage tanks.

VERSION I




VERSION II


The best results in „Low Volume” Tank designs

All Colenta Processors are equipped with precision temperature control of all processing chemistry with automatic Developer Tank cooling to prevent overheating. Replenishment is preheated by injection into the inline chemical heating systems, which are controlled to pre-set programmed temperatures within +/- 0,2 °C.

TRANSPORT SYSTEM:

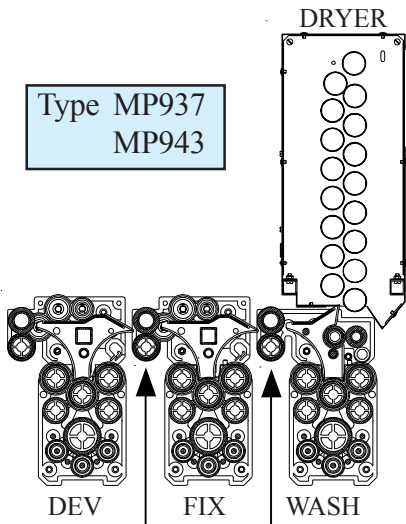

Rinsed cross over roller

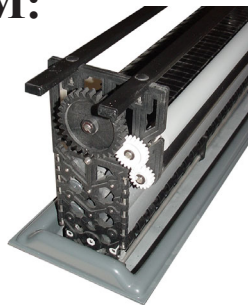


Optional accessories:

Rack Carrier Tray/Drip Tray for easy handling the transport racks.

Cover Lid for Film Feed Tray



Closed stand for processor

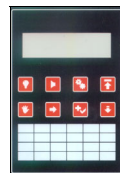

Microprocessor controlled Displa



Technical specification subject to change without notice.

see additional Data sheet: N1-15

COLENTA LABORTECHNIK
Ges.m.b.H. & Co. KG
Neunkirchner Strasse 117
A-2700 Wiener Neustadt
Tel.: +43-2622-28311-0
Fax: +43-2622-28311-7
E-mail: office@colenta.at
<http://www.colenta.at>

ACI-COLENTA
INTERNATIONAL Div.
Forge Close, Little End Rd,
Eaton Socon, Cambs. UK
Tel. +44-1480-211118
Fax: +44-1480-215236
Fax: +44/870/132 1352
E-mail: aci@patflaherty.com