



# Colenta<sup>®</sup> MEDIPHOT 903

The **Colenta MP903** is designed as a free standing processor for highest processing needs. Made for large hospital/clinics operating at 24/7. The **MP903** offers a capacity of up to 300 films/hour size 35x43 cm film at a 90 sec. Dry-to-Dry fully automatic & monitored processing cycle.

## Technical Specifications MP902:


max. Film width:	43cm/17inch
min. Film length:	10cm/4inch
Tankbody PVC:	Dev/Fixer/Wash 10 lt. each
Dev. Temp. Range:	18-43°C ± 0,2°C
Dryer Temperature:	monitored at 50°C
Replenishment:	film based / programmable for D&F
processing capacity	@ 90 sec Dry to Dry time: <b>175cm/min</b>



The **Colenta MP903**, as all Colenta medical X-Ray processor models, provides a fully automatic „dry to dry“ film processing technology. A microprocessor system ensures precise temperature control of all processing chemistries with automatic Developer / Fixer tank cooling to prevent overheating. Replenishment is preheated by injection into the inline chemical heating systems, which are controlled to pre-set programmed temperatures within +/- 0,2 °C.

## Standard Features and accessories:

- Low Level monitoring of all Tank solutions.
- Automatic Solution Tank „Top-Up“ control on „Start-Up“
- „Work in Progress“ monitoring
- inbuilt repl. storage tanks (2x30L)
- prepared to connect to be connected an external air exhaust system
- drive motor control
- water rinses cross over rollers
- small spare part kit
- CE marking

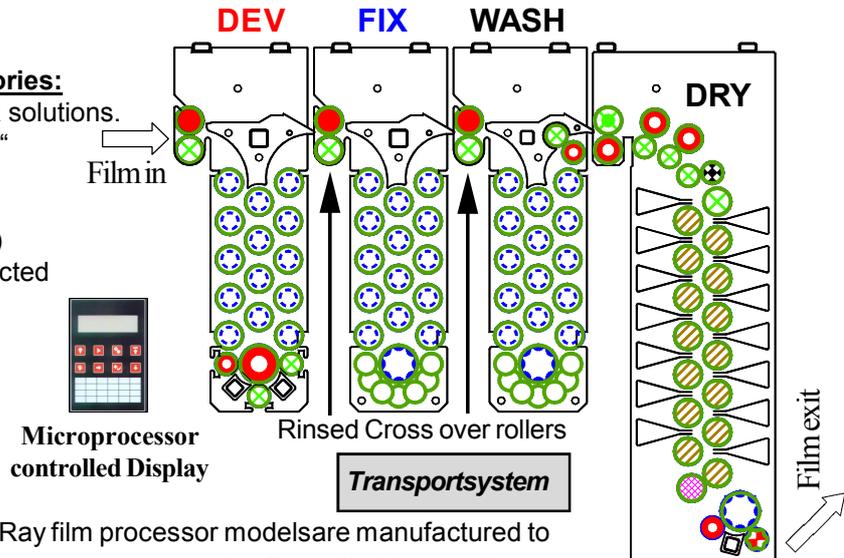


All **Colenta Mediphot** medical X-Ray film processor models are manufactured to meet the highest standards of health and safety, are all CE certified in compliance to the strictly controlled European laws that are applied to the production and operation of our products.



## optional accessories:

- 1) Cover Lid for Film Feed Tray
- 2) chemical mixing console...see extra PDF
- 3) water supply panel..see extra PDF
- 4) rack drip tray

.....  
more see Mediphot Price list!

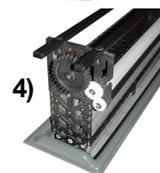

# Colenta<sup>®</sup>

**COLENTALABORTECHNIK**  
Ges.m.b.H. & Co. KG  
Neunkirchner Strasse 117  
A-2700 Wiener Neustadt  
Tel.: +43-2622-28311-0  
Fax: +43-2622-28311-7  
E-mail: [office@colenta.at](mailto:office@colenta.at)  
[www.colenta.at](http://www.colenta.at)

## basic installation requirements MP903:

Power Supply:	3 x 230/400 VAC + N , 50/60 Hz, 5,5kW
Water Supply:	R3/4", max 15°C, 3-6 bar
Water consumption	2,5 litre/min when film in proc
Water Drain:	32mm (5/4inch) hose or tube

**See additional, technical data sheet & preinstallation manual.**